The background features a repeating pattern of overlapping circles and stylized leaves. The circles are arranged in a hexagonal grid, with each circle overlapping its neighbors. Inside each circle, there are two stylized leaf shapes, one on the left and one on the right, pointing towards the center. The entire pattern is rendered in a light gray line art style on a white background.

LED FACE MASK

the beauty specialists

ENGLISH	3	ESPAÑOL	31
FRANÇAIS	10	PORTUGUÊS	38
DEUTSCH	17	ITALIANO	45
NEDERLANDS	24	УКРАЇНСЬКА	52

WARNINGS & PRECAUTIONS

- **Please read all instructions thoroughly and retain for future reference.**
- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- If the power lead is damaged, it must be replaced by the manufacturer or its service agent or similarly qualified person in order to avoid hazard. Contact Customer Care to arrange a replacement.
- Do not immerse product in water. This appliance may be cleaned by wiping with a damp cloth. Do NOT use solvents or abrasive cleaners.
- Do not use the LED face mask if you are taking any medication that may cause sensitivity to light.
- Using the mask while on light-sensitizing medication may lead to skin irritation or other reactions.
- Some medications, such as certain antibiotics, cancer treatments, and diuretics, can make your skin more sensitive to light, so always check with your doctor if you are unsure.
- Use the device only for the specific skin conditions listed in the instructions and not for any other purpose.
- Do not use the mask if you are pregnant, breastfeeding, or trying to become pregnant.
- Avoid using the device if you have a medical condition like Lupus, Albinism, or photosensitive eczema.
- Do not use the mask if you have any disorder that causes sensitivity to light.
- Always perform a patch test before using the device for the first time to check for any skin sensitivity.
- DO NOT use LED Face Mask in any way that is not exactly as specified in these instructions for use.
- Stop using the mask immediately if it feels hot or causes your skin to feel uncomfortably warm.

- Use only the recommended gel masks that are designed to work safely with this device.
- Make sure to clean your face and remove any makeup, creams, or lotions before using the mask.
- Keep the device away from water and do not use it in wet environments or if it may have come into contact with liquids.
- Do not place the device near hot surfaces or expose it to high temperatures.
- Never try to modify or repair the device yourself as doing so can void the warranty and may be unsafe.
- Do not use the mask if any part such as the mask, controller, or charging cable is damaged.
- The strap contains polyester, nylon, and lycra, so do not use it if you are allergic to any of these materials.
- Do not use the device if you are prone to headaches triggered by bright or flashing lights.
- Always attach the correct adapter before plugging the device into a socket
- DO NOT share LED Face Mask with any other person. LED Face Mask is for single person use only.
- Follow the recommended treatment guidelines and avoid using the mask more often than advised.
- LED Face Mask must only be supplied at safety extra low voltage (SELV) corresponding to the marking on LED Face Mask.
- Do not use on areas that have had cosmetic surgery.
- Always unplug the device immediately after using and before cleaning.
- Do not leave the device unattended. Turn off before you leave the room.
- Keep out of reach from children and pets.
- Do not use in the bath or shower.
- Do not charge for longer than 24 hours.
- LED Face Mask is only to be used with the controller supplied (LED Face Mask Controller).
- If you suffer from a skin condition, consult your dermatologist before use;

Sensitivity to light, or photosensitivity, is a known side effect of many medications. This may include certain types of antibiotics, cancer treatments, and water pills (diuretics). If you're not sure whether a medication you're taking could cause this reaction, it's best to check with your doctor or pharmacist.

In addition to medications, other substances can also make the skin more reactive to sunlight. Some common examples include: St. John's Wort, coal tar, some deodorants, antibacterial soaps, artificial sweeteners, mothballs (which contain naphthalene), petroleum-based products, optical brighteners in laundry detergents, and cadmium sulfide (sometimes used in tattoo pigments).

INTRODUCTION

With natural aging, the skin gradually loses its firmness, moisture, and elasticity, often leading to wrinkles and dullness. This LED Face Mask is designed to support facial rejuvenation using targeted wavelengths of light therapy.

The mask emits light at specific wavelengths which are known to support the reduction of fine lines and wrinkles and promote healthier-looking skin over time.

This device is intended for home use by adults seeking a non-invasive option to improve facial appearance and reduce the visibility of fine lines. Results may vary based on individual skin condition and usage frequency.

KNOW YOUR UNIT

- | | | |
|----------------------|-------------------------|---------------------|
| ❶ Head Straps | ❸ Mode/Colour Indicator | ❾ Battery Indicator |
| ❷ LED Face Mask | ❹ Intensity Indicator | ❿ Time |
| ❸ USB Charging Cable | ❺ On/Off | ⓫ Mode |
| ❹ Carry bag | ❻ Intensity Increase | ⓬ Time Increase |
| ❺ Controller | ❼ Intensity Decrease | ⓭ Time Decrease |

PATCH TEST

Before beginning your first full-face treatment, it's important to conduct a sensitivity test:

1. Choose a clean, small area of skin (such as your inner forearm).
2. Place the LED mask over this area and run a complete session at the planned treatment intensity level.
3. Wait 6–8 hours after the test. If you notice excessive redness, itching, or irritation, do not proceed with facial use.
4. If no adverse effects are observed, you may begin normal treatments.

CHARGING


To charge the LED Face Mask, insert the USB Charging Cable into the Controller and plug into a USB port. The bars on the Battery Indicator will show 5 bars once fully charged. The Controller will take between 2.5-3 hours to fully charge.

Note: Always recharge when the Battery Indicator has 1 bar. The LED Face Mask will not operate whilst charging.

PREPARING FOR USE

Before each session, ensure the LED Face Mask is powered off and clean. Refer to the cleaning section on the next page for detailed instructions. The mask is equipped with adjustable head straps for a secure and comfortable fit:

1. Insert the straps through the slots located on each side of the mask.
2. Wrap the straps around your head and fasten them using the provided hook-and-loop system. Adjust for comfort and stability.
3. Insert Lead from the LED Face Mask into the Controller.

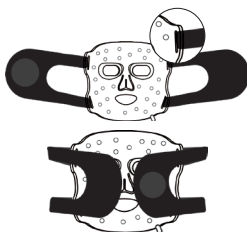TREATMENTS

To achieve optimal results, it is recommended to use the LED Face Mask consistently over several weeks.

- Begin with 3 to 5 sessions per week for the first 6 weeks.
 - Adjust frequency based on your skin's sensitivity and your skincare routine.
1. Press the **On/Off** button to turn on the mask. The default mode (**Mode 1**: Red + Near-Infrared Light) will activate automatically.
 2. To select a light mode, press the **"Mode"** button to cycle through the four treatment options:
 - Mode 1**: Red + Near-Infrared Light (17minutes maximum)
 - Mode 2**: Red Light Only (25minutes maximum)
 - Mode 3**: Blue Light (20minutes maximum)
 - Mode 4**: Red + Blue Light (20minutes maximum)
 3. To adjust the intensity, use the **"Intensity +"** or **"Intensity -"** buttons to set brightness (Levels 1–3). Start at a lower intensity and increase based on comfort.
 4. To set a treatment time, press the **"Time +"** or **"Time -"** buttons to adjust session duration in **5-minute increments**. The timer screen will countdown and the mask will automatically turn off once the treatment time has finished.

Note: Each mode has specific maximum time limits - see above.
 5. Once the session ends, remove the mask, turn off the device, and apply your regular moisturiser and sunscreen

AFTERCARE & REACTIONS

Following a treatment, some individuals may notice temporary skin responses:

- Slight redness can appear post-use but typically fades within a few hours.
- In some cases, users may experience minor peeling, dryness, or itchiness, which usually resolves in a short time.
- Occasional mild headaches have been reported after sessions.

The device uses light wavelengths generally considered safe for cosmetic use. However, if you experience any unusual discomfort or skin reaction, discontinue use immediately and consult a healthcare professional if necessary.

Optional Use of Performance Enhancing Gels:

The benefits of using the LED mask can be boosted if used in conjunction with specific proprietary approved gels. For more information visit www.riobeauty.com.

CLEANING

Proper cleaning is essential to maintain hygiene and prolong device life.

Before each use:

- Power off the device and unplug from any charging source.
- Wipe the interior and exterior of the mask with a soft, damp cloth.
- Avoid harsh chemicals, abrasive materials, or submerging the unit in water.
- For deeper cleaning, lightly use a cloth dampened with 70% alcohol.

Controller Care:

- Wipe gently with a dry or slightly damp cloth.
- Avoid moisture exposure near connection ports.

STORAGE

- Store the mask in a clean, dry place at room temperature.
- Avoid exposure to direct sunlight, excessive humidity, or dust.
- Use the provided carry bag to protect the device during travel or when not in use.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No lights turn on	Battery not charged	Fully charge the device using the supplied USB cable.
Lights stay on beyond timer	Timer function error	Turn off manually. If this repeats, contact customer support.
Flickering lights	Internal malfunction	Discontinue use and request a replacement.
Device won't turn on	Loose connection	Check all cables and ensure a full charge. Retry after disconnecting.

Note: If a problem persists after troubleshooting, stop using the device and contact Rio Customer Care.

DISPOSAL

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

Recycle Responsibly:

1. Check your local municipal guidelines
2. Separate the components and dispose of them correctly

SPECIFICATION

Controller Input: 5V  1.0A 5.0W 

Mask Input: 5V  1.0A 5.0W Max from controller

This product complies with current EU product legislation.

This symbol indicates the USB connector.

Read these instructions carefully before use and then retain for future use.

Class III symbol (IEC 60417-5180).

This product is for indoor use only.

This unit complies with current UKCA product legislation.

Direct current.

Recycle Paper.

Detachable PSU IEC 60417-6181.

Model number: FCAA

Program time: 5-25 minutes

Wavelengths: • Red (650-655nm)
• Blue (405-410nm)
• NIR (845-855nm)

Maximum output intensity:

30mW/cm²

Operating temperature:

0°C to 30°C (32°F to 86°F)

Storage and transport temperature:

-5°C to 55°C (23°F to 131°F)

Operating humidity:

30% to 95% without condensation

Storage and transport humidity:

20% to 85% without condensation

Electrical output: 5V, 1A

IP22

CUSTOMER CARE

If you have any questions, please contact us at: online@riobeauty.com

To reorder additional accessories, please visit
www.riobeauty.com or contact your local distributor.

GUARANTEE

To register your two year guarantee online, visit:

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©The Dezac Group Ltd 2025, PO Box 17, Cheltenham Spa, England GL53 7ET
UK RP: The Dezac Group Ltd, PO Box 17,
EU RP: The Dezac Group (Europe) B.V. Azura,
UK Product RP: The Dezac Group Ltd, PO Box 17, Cheltenham Spa, England, GL53 7ET
EU Product RP: The Dezac Group (Europe) B.V. Azura, Saturnusstraat 46-62, 2132 HB Hoofddorp, Netherlands